

**TANITA**

Monitoring Your Health

# *FacePeek!*<sup>TM</sup>

Instruction Manual  
*Manuel d'instructions*

## Introduction

Thank you for purchasing FacePeek™ Skin Condition Monitor. For maximum benefit, please take a few moments to read through this manual, and keep it on hand for reference. With FacePeek, you're on your way to beautiful skin.


### *Where Beautiful Skin Begins*

FacePeek uses a harmless, non-invasive, handheld sensor system to monitor the environment at the surface and below where new skin starts forming. It provides a relative measure of:

- softness
- moisture
- oiliness

FacePeek sends out safe ultrasonic waves to measure the relative condition of your skin, from the surface to the tissue approximately 1.5 mm below. Frequent monitoring along with appropriate care (see your enclosed guide: *Skin Care at a Glance*) can be helpful in maintaining or improving your complexion and appearance.


### *FacePeek Sensor System*

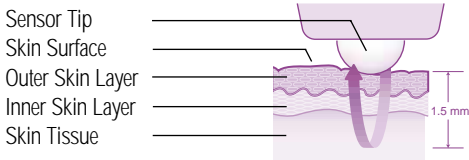

# System Components/Functions



- |               |                        |
|---------------|------------------------|
| 1. Cap        | 5. Advance Button      |
| 2. Sensor Tip | 6. On/Off Button       |
| 3. LCD Panel  | 7. Battery Compartment |
| 4. Set Button |                        |

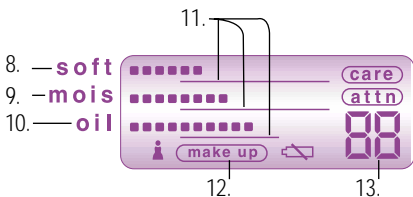


## LCD Panel Elements

- |                             |                         |
|-----------------------------|-------------------------|
| 8. Relative Softness Level  | 11. Typical Level Lines |
| 9. Relative Moisture Level  | 12. Makeup Setting Icon |
| 10. Relative Oiliness Level | 13. Age Setting         |





# Operating FacePeek
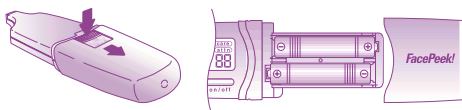
## *Inserting Batteries*

1. Remove the cap and battery cover.
2. Insert the batteries. Take care to insert the batteries to the correct polarity.
3. Slide on the battery cover and replace cap.



## *Inputting Data*

1. Press the On/Off power button.
2. Wait five seconds for all indications to appear on display and unit to beep.
3. Press the Set button.
4. LCD panel will flash the age icon.
5. Press the Advance (▲) button to input your age.
6. Numbers will increase each time you press.
7. To fast forward, hold the button down. (After age 70, system will recount from 0.)
8. Press the Set button to record data.
9. Press the Advance (▲) button to select makeup icon:
  - female without makeup 
  - female with makeup  **make up**
  - male without makeup 
  - male with makeup  **make up**
- (Makeup is defined as cosmetics, moisturizer, or sunscreen.)
10. Press the Set button to record data.

## Taking Measurements

1. Place the sensor tip under your cheekbone.
2. Position the unit at a right angle (90°) and press gently to skin.
3. Hold there until you hear a beep.
4. Measurement has now been completed.

## Reading Results

Your relative measurements are indicated by dots displayed above the Typical Level lines (shown in purple) which represent customary ranges.

**s o f t**    ●●●●●●●●●●●●●●●●●●●●

**m o i s**    ●●●●●●●●●●●●●●●●●●●●

**o i l**        ●●●●●●●●●●●●●●●●●●●●

**Soft**    Dots indicate degree of softness, from low to high.

**Mois**    Dots indicate degree of moisture, from low to high.

**Oil**      Dots indicate degree of oiliness, from low to high.

If your results fall within the length of the Typical Level lines, your skin is within customary range. If your results fall below or beyond, special icons appear:

**care**    appears when your relative results are not within typical range.

**attn**    appears when your relative results are unusually high or low.

## To Obtain Optimal Results

- Position the sensor at a right angle under your cheekbone: the Typical Level lines are based on measuring in this area. However, it is possible to obtain readings in other areas, such as your "T-Zone" (forehead and nose).
- Apply minimal pressure as you gently touch sensor to skin.
- If your skin is sticky, remove makeup and try again.
- If your skin is dry, apply foundation and select "makeup" mode.
- See handling tips to properly care for unit.

## About Your Skin

### *Softness*

Levels of softness say a lot about inner skin balance. Moisture and oiliness are major components. Softness eases tightness and can help minimize lines and wrinkles.

### *Moisture*

Levels of moisture are related to changes in the sebaceous membrane protecting keratinous cells. Moisture contributes to softness, and can alleviate dryness and irritation. If levels are low, skin is less elastic.

### *Oiliness*

Levels of oiliness are related to secretions of sebaceous matter in the membrane protecting keratinous cells. A normal amount is needed to maintain moisture. Excess oiliness can cause complexion problems. Most commonly, oily skin appears in the "T-Zone".

Many factors such as climate, diet, medication, stress, and cosmetics can affect skin condition. For steps to healthier skin, please refer to your enclosed guide from *FacePeek: Skin Care at a Glance*.

## More About FacePeek

### *Handling Tips*

- Do not disassemble or wash the unit.
- Frequently clean off cosmetic residue with dry cotton swab.
- Do not use water, benzene, thinner, alcohol or other liquids.
- Store with cap on, free from heat, humidity, direct sunlight, and dust.
- Do not apply heavy pressure or vibration to the unit.

### *Troubleshooting*

**If nothing appears on the display after pushing the power button:**

- Check to see if batteries are correctly positioned.
- If batteries are low, replace them.

**If readings vary every time you measure:**

- Place the sensor tip to the same area under your cheekbone every time.
- Hold the unit at a right angle every time.
- Apply a little more pressure.
- Clean the sensor with a dry cotton swab.

## Replacing Batteries

The LCD panel will flash a Battery Icon to indicate when batteries are low. When this occurs, replace both batteries (AAA type) together.

FacePeek comes with batteries; their service life may be short. Replacement batteries can last for about one year if the unit is used approximately once a day.

If you do not intend to use FacePeek for an extended period of time, remove batteries to prolong life.

### Notes

- Individual skin care needs will vary. Consult your physician or dermatologist about your specific skin type and your skin care options.
- Discontinue use of any skin care product if rash develops or if you experience discomfort. Consult your physician or dermatologist if discomfort continues.
- Information provided as a guide only. Results may vary depending on skin type and condition and Tanita does not guarantee or warrant improvement in appearance in skin condition or complexion.

# Warranty

## Not Valid Without Proof of Purchase

Your Tanita product is warranted to the original consumer only, for one (1) year, to be free from defects in materials and workmanship under normal use. This product is designed **for personal use only**. Tanita Corporation of America reserves the right to nullify the warranty agreement should the product be used in a commercial environment. This warranty is void if the product has been tampered with or mistreated, or if the defect is related to servicing or parts not supplied by Tanita. Tanita's obligation under this warranty is limited solely to the charge of repairing the product, plus the cost of parts or, at the option of Tanita, replacing the product with an equivalent model. There are no express or implied warranties except as listed above. **Tanita is not liable for incidental or consequential damages for failure or malfunction of the product and any recovery may not be greater than the purchase price paid for this product.** Some states do not allow limitation on how long an implied warranty lasts, or the exclusion or limitation of incidental or consequential damages. This warranty gives you specific legal rights. These rights may vary from state to state. **No person is authorized to change the terms of this warranty.**

## If Product Is Defective or Needs Repair

Do not return this product to the place of purchase; doing so may delay the processing of your warranty claim.

Call Tanita at 1-800-TANITA-8 to receive a Return Authorization Number before returning this product. **Proof of purchase must accompany the product.** Because of possible loss, we recommend insuring your product, return receipt requested, when using the mail. If you do not receive a receipt within a reasonable time, start a tracer through the originating post office. It is important that you pack the product carefully prior to shipment to prevent damage in transit and that you include a letter describing in detail the specific condition requiring adjustment or repair.

## After Warranty Repairs

Please call Tanita at 1-800-TANITA-8 (USA and Canada only).

Please complete the following information and retain for your records.

**Model: FacePeek HC-220**

---

**Date Purchased:**

---

**Purchased From:**

---

**Address:**

---

---

# Specifications

## Model

FacePeek (HC-220)

---

## Display Results

3 levels

---

## Dimensions

5" H x 1-1/2" W x 7/8" D

---

## Battery

2 x AAA (included)

---

Should you need to return this product for an exchange or refund, please contact the store/cataloger from which it was purchased. Or call 1-800-TANITA-8 (USA and Canada only).

# Spécifications

## Modèle

FacePeek (HC-220)

---

## Résultats d'affichage

3 niveaux

---

## Dimensions

5 po H x 1-1/2 po L x 7/8 po P

---

## Piles

2 x AAA (incluses)

---

Si vous devez retourner le présent produit pour un échange ou un remboursement, veuillez communiquer avec le magasin/catalogue où vous avez acheté le produit ou téléphonez au 1 800 TANITA 8 (États-Unis et Canada seulement).

---

# TANITA

Monitoring Your Health

[www.tanita.com](http://www.tanita.com)

Tanita Corporation of America Inc.  
2625 South Clearbrook Drive  
Arlington Heights, IL 60005, USA  
Tel: 847-640-9241  
Fax: 847-640-9261

© 2001 Tanita Corporation.  
All rights reserved