

**MOBI®**

User Manual v1.3  
Item 70678

# Smart Medication & Vitamin **Wi-Fi Pill Dispenser** With UV Blocking Lid



**User Manual**

Thank you for selecting the MOBI® Smart Medication & Vitamin Wi-Fi Pill Dispenser. Our app simplifies the process of connecting your device for convenient health monitoring. You can effortlessly transmit your measurements wirelessly to your mobile device and access them through the MOBI® Smart App.

## Table Of Contents

Table of Contents .....	2
What's In the Box .....	3
Product Information .....	3
Operating Device(Without the App) .....	6
Operating Device(With the App) .....	7
Features .....	24
Specifications .....	26
Care & Maintenance .....	30
FCC Statement .....	31
Limited Product Warranty .....	34
Extended Manual .....	35




## What's In the Box

1. 1 x Pill Dispenser
2. 9 x Dose Rings
3. 1 x Locking Key
4. 1 x Power Cable
5. 1 x Power Adapter
6. 4 x AA Batteries
7. 1 x User Manual

## Product Information

Stay at ease with the MOBI® Smart Medication & Vitamin Wi-Fi Pill Dispenser. Through the MOBI® Smart App, you can set up to 6 alarms per day and have 4 different alert systems to help you keep track of the dosages that must be taken. With its automatic detection and rotation, you can be sure that you are taking the correct dosage at the correct time.

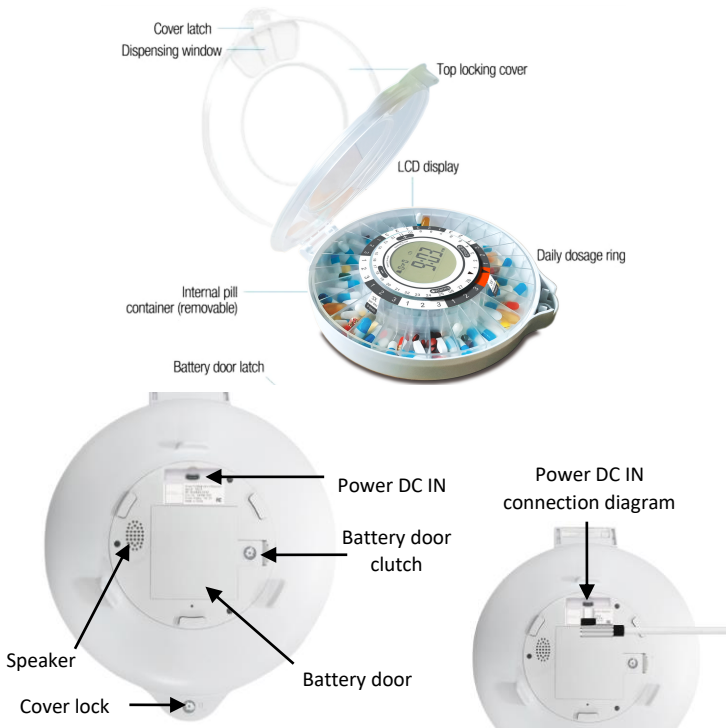

Your Smart Pill Medication Dispenser comes with **9 Daily Dosage Rings** corresponding to how frequently medication needs to be taken daily.



## Device Setup

1. Use the supplied Security Key to install the provided backup batteries.
2. To access the LCD Display touch screen, open the top cover.
3. Set the Hour and Time by pressing the "CLOCK" button.
  - a. Press the "+" and "-" buttons to adjust the hours & minutes. Press "CLOCK" each time to save.
4. Set alarms by pressing the "ALARMS" button and using the "+" and "-" buttons to choose an alarm.
  - a. Press "ALARMS" each time to save/exit. Press the "+" and "-" buttons to adjust hours & minutes.
5. Set the Alarm Tone by holding the "ALARMS" button for 3 seconds to choose an Alarm Tone.
  - a. Hold the "-" button for 3 seconds to set the Alarm Volume level (silent, low, medium, & high)
6. Press & hold "CLOCK" & "ALARMS" buttons for 5 seconds to reset it to factory settings.

## Operating Device (Without the App)

1. Place the internal pill container inside Smart Pill Medication Dispenser.
2. Choose the daily dosage ring that corresponds to your needs.
3. Place the daily dosage ring in the center of the internal pill container and fill your medications in the slots corresponding to the daily dosage ring.
4. Set the clock, alarm volume, alarm type, and alarm times.
5. Close the cover and lock using the security key. (Please keep key in a safe & accessible location)
6. Close the medication dispenser window and place Pill Dispenser on a flat surface, out of direct sunlight.
7. To dispense medication and to turn off the Alarm / Notification
  - a. Open the small dispensing window on top.
  - b. Lift & tilt the pill box dispenser.
  - c. Slide to close dispensing window.



## Operating Device (With the App) Download the MOBI® Smart App



- Ensure you have access to your email to verify your MOBI® Smart account.
- Have your 2.4 GHz Wi-Fi Network Name & Password ready.
- For the best experience, please grant the MOBI® Smart App access to the following Permissions:
  1. Location Access
  2. Local Network Access
  3. Bluetooth Access



## Register

Open MOBI® Smart App and tap Sign Up. Agree to the privacy policy. Enter your email address, then tap “Get Verification Code”. Check your email for the verification code and then create your password.



## Add Device

Tap “Add Device”, then select your MOBI® Smart Device. (You may also tap the plus sign ‘+’ on upper right-hand corner to add a device).



App screen will auto show the device able to be Chosen.



## Enter Wi-Fi

"Connect your mobile device to a 2.4GHz Wi-Fi network via the phone's settings. Then within the app, enter/select the same 2.4 GHz Wi-Fi network and tap next.

**Note:** Wi-Fi name and password are case, space and number sensitive".

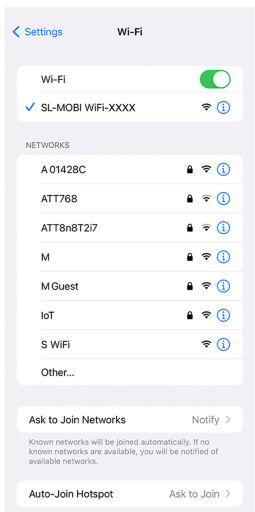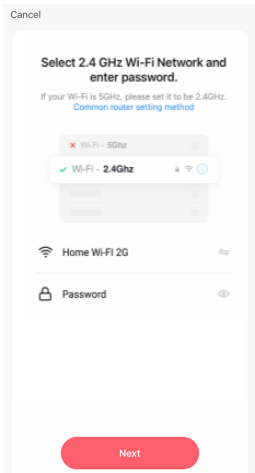


## Connect to Hotspot

Follow the in-app instructions to reset the pill dispenser. Once the unit is reset, click "Connect" to be redirected to your Wi-Fi settings. Within the Wi-Fi settings please select the available network appearing as "SL-MOBI Wi-Fi-XXXX". After connecting, return to the MOBI® Smart app to continue.

## Success

Added.

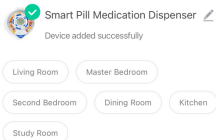

If your device has been successfully paired, you'll hear a voice confirmation from the Smart Pill Medication Dispenser. Once connected, you may rename your device for easier management, or edit later by clicking the pencil icon. Tap "Done" to be taken to the Hub controls.

Scan this QR code to register your product on our website for a FREE warranty extension.

If you experience any issues, please contact us at [www.mobiusa.com/gethelp](http://www.mobiusa.com/gethelp).




You can also reach us at +1 (323) 904-3000 between Monday through Friday, 8:30 AM - 5:00 PM PST.



## Understanding of Daily Dosage Rings

Your smart pill dispenser comes with 9 Daily Dosage Rings corresponding to how frequently medication needs to be taken daily.

The rings help properly dispense medication from the internal pill container and include 28 daily medication dosage slots.

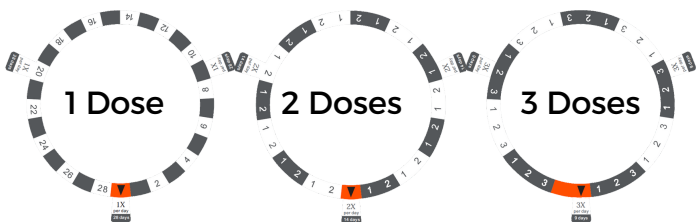


First, place the dosage ring in the center of the internal pill container with the tab facing the bottom (indicated by the ▼ icon on the plastic container.)

Then align the two triangle marks on the dosage ring and faceplate, then fold the tabs to secure the dosage ring before filling the dispenser with medication.

Once the dosage ring is correctly situated, the alternating colors help distinguish between different days. For example, if you use the 2x per day ring (meaning medication only needs to be dispensed twice daily), the internal pill container will hold 14 days' worth of medication.

The grey spaces to the right of the orange tab represent 1 days' worth of medication. The following white spaces represent the upcoming 2<sup>nd</sup> day's medication. The dosage rings will alternate between grey and white spaces to distinguish between different days. Each colored space includes numbered slots representing the different alarm times configured during the device setup process.

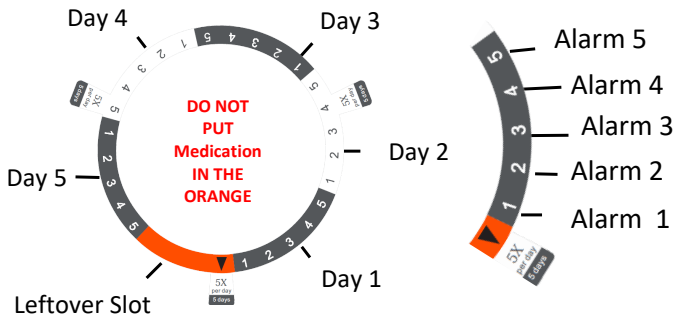

If you need medication dispensed 5 times during the day, the ring can accommodate 5 days or 25 doses, with three leftover slots.

Leftover slots are indicated by the orange spaces on the ring and are substantially larger orange areas on the 5x per day and 6x per day rings.

**Reminder: Do not put medication in leftover slots.**

Choose the Daily Dosage Ring with the tab that shows how often you need medication dispensed daily. Position the ring around the compartment. Fill each compartment accordingly, then place the carousel back into the dispenser with the dosage ring.

## Plugging in the Cable

1. Connect the power cable to the pill dispenser.
2. Plug the other end into a USB port.

**NOTE:** The cable should always be plugged in, as it is the primary power source for this device.

Ring	# of Meds per Day	Programable Table Period
1	1	28 Days
2	2	14 Days
3	3	9 Days
4	4	7 Days
5	5	5 Days
6	6	4 Days
7	7	4 Days
8	8	3 Days
9	9	3 Days



## Installing the Batteries

1. Turn the dispenser upside-down. Insert the locking key into the battery cover and turn clockwise.
2. Remove the battery door.
3. Insert four (4) AA alkaline batteries inside the compartment.
4. Replace the battery cover. Insert the key into the cover and turn counterclockwise to lock.


**NOTE:** Replace batteries if the low battery symbol appears.

Batteries are installed as backup power in case of power outages. The device will stay on for approximately (4) days when powered by batteries.

## Locking and Unlocking the Cover

The cover lock helps prevent accidental spills and unauthorized access and can also help with confused patients who may try to take medication incorrectly.

1. Align the locking key with the lock hole (as shown.)
2. **To lock;** turn the key clockwise while pressing down on the cover latch.
3. **To unlock;** Turn the key counterclockwise.

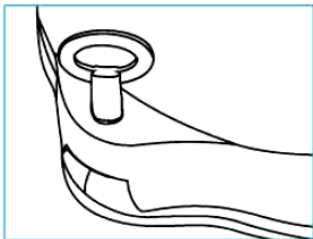

Lift the cover open using the cover latch.


## Selecting 12-Hour/24-Hour Mode

Press and hold the CLOCK button for (2) seconds to switch between 12 and 24-hour mode.



## Setting the Time

1. Press the CLOCK button once to start the time setting.
2. Press the + or - buttons to adjust the hour. Then, press the CLOCK button again to adjust the minute.

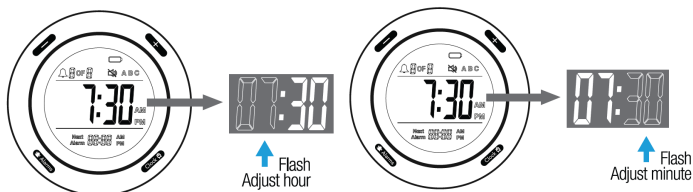


3. Use the + or - button to adjust the minute, similar to the hour.

**NOTE:** Pressing and holding the + or - button will cause the hour or minute digits to move faster.

4. Press the CLOCK button to exit time settings with changes saved or press and hold the CLOCK button for two (2) seconds to cancel the setting change and exit out of the time setting.

**NOTE:** The digits on display will flash to indicate which digits are currently being changed.


## Setting the Alarm

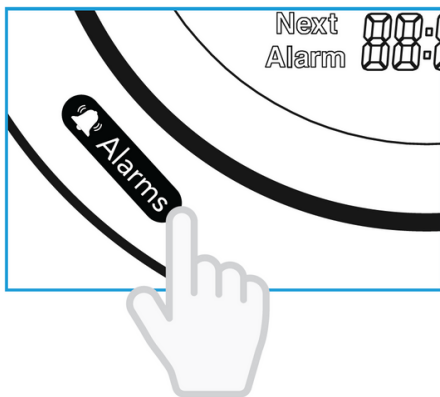

1. Press the "ALARM" button to start the alarm time setup.
2. Press the "+" or "-" to select between an alarm, then press the "ALARM" button to confirm.
3. Press the "+" or "-" to adjust the hour. Then, press the "ALARM" button to confirm the hour and set the minute.
4. Press the "+" or "-" button to adjust the minute, similar to the hour.




**NOTE:** Pressing and holding the "+" or "-" button will cause the hour, minute digits to move faster. The digits on display will flash to indicate they are currently being changed.

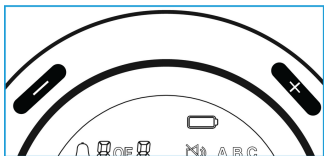

5. Press the "ALARM" button once or press and hold down the "ALARM" button for two (2) seconds to save and exit alarm setting mode.

To cancel the setting change, wait a few seconds for alarm setting mode to end automatically.

To delete current alarm settings, press and hold the "CLOCK" button while in alarm setting mode until  is shown on the display screen.

## Setting the Alarm Tone & Volume

1. Press and hold the "ALARM" button for three (3) seconds to select alarm tones. On-screen alarm tones are displayed as A, B, and C.
2. Press and hold the "-" button for three (3) seconds to change volume levels.


**NOTE:** Volume level is represented by On-Screen "speaker" symbol.

The three (3) alarm tones have distinct beeping patterns:

- Alarm A: 20 seconds of continuous beeping, followed by 10 seconds of silence, then 20 seconds of continuous beeping as the pattern repeats.
- Alarm B: Three sets of three beeps. The beeps will stop for 7 seconds, and then resume beeping, as the pattern repeats.
- Alarm C: One (1) beep every 3 seconds.




**NOTE:** All alarms will continue their patterned beeping for 30 minutes before stopping. Also, all alarms will sound with the same alarm, until it is manually changed to a different tone.

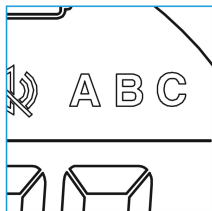

## Reset Factory Settings

Press and hold the "CLOCK" and "ALARM" buttons simultaneously for five (5) seconds to clear and reset all settings.

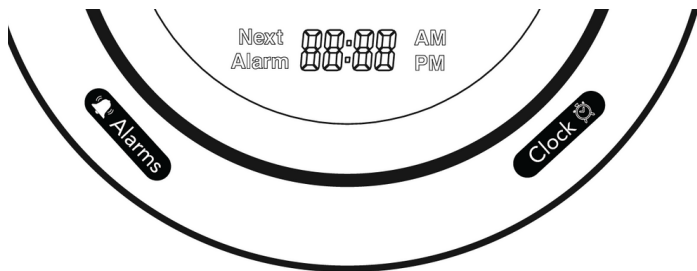

## Operating the Smart Medication & Vitamin Wi-Fi Pill Dispenser

1. Choose the Daily Dosage Ring that corresponds to your needs (please refer to "Understanding the Daily Dosage Rings" on page 11).
2. Situate the Daily Dosage Ring in the center of the internal pill container and place your medications in the container to correspond with when they are taken (1st, 2nd, 3rd, etc.) indicated by the numbers on the ring (1, 2, 3, etc.) Do not place any medication in the orange sections. When the orange is approaching the dosage window, it is time to refill the container.
3. Install the batteries (refer to "Installing the Batteries" on page 15).
4. Unlock the cover using the key (please refer to "Locking and Unlocking the Cover" on page 15) and place the internal pill container inside the main body.
5. Set the time (please refer to "Setting the Time" on page 17).



6. Adjust the alarm volume and alarm tones and choose a 12-hour or 24-hour time (please refer to "Setting the Alarm Tone and Volume" on page 20) and ("Selecting the 12-Hour/24-Hour Mode" on page 16.)
7. Set your medication times (refer to "Setting the Alarm" on page 18).
8. Close the cover and lock it with the provided key (please refer to "Locking and Unlocking the Cover" on page 15).
9. At designated medication times, the alert will sound.
10. Slide open the window by pushing the tab to the right.
11. Tilt the pill dispenser to pour the pills into your hand or medication cup. Make sure all pills have come out of the compartment. The alarm should stop automatically.
12. Close the dispensing window and place the pill dispenser on a flat surface.

## Features

- Using the MOBI® App, you can set up to 6 alarms per day to help you keep track of daily medicine dosage that has to be taken. Receive app notifications and text alerts.
- The Pill Dispenser also has a Flashing LED light and a Voice Monitor that alerts you when it is time to take your daily dosage.
- The MOBI® Pill Dispenser has a key to physically lock the device, so you don't have to worry about medicine abuse.
- The tray is replaceable and easy to remove for cleaning.
- It has 28 slots with the capacity to hold 35 pcs of Aspirin pills.
- Large text for easy reading.
- Automatic detection to see if medication has been taken.
- Low Battery Indicator.
- Compatible with IOS and Android devices.
- Automatic rotation when it is time to take medicine.
- Device settings can be set manually or through the MOBI® app.



- Wi-Fi connection allows for multiple family members near or far to share in monitoring for the patient.
- Works with the MOBI® Smart App along with other MOBI® Health devices to provide more comprehensive and managed caregiving to your loved ones.
- Create Automations with other MOBI® Smart products as part of the Home Health Environment.
- Allows patients & caregivers to sort, schedule, dispense, monitor, and share data. Automatically dispense medications for up to 28 days with tracking, smart alerts, and notifications. MOBI® smart app helps monitor your health with other MOBI® devices all in one app!


## Specifications

1. Dispenser length: 9.25"
2. Dispenser width: 9.25"
3. Dispenser height: 3"
4. Weight: 1.9lbs.
5. Power: 4x AA alkaline batteries (included)
6. Product type: Pill dispenser
7. Approved Location: Indoor use only
8. Adult assembly required: No
9. Additional tools required: No
10. Number of compartments: 28
11. Number of dosage rings: 9
12. Connectivity: 2.4 GHz Wi-Fi
13. Alarms: Audible and visual
14. Cover lock: Yes
15. Country of origin: China
16. Warranty: Yes



## Safety Tips & Precautions

- This instruction manual is meant to ensure the safe operation of the product. Non-compliant use may void the warranty.
- Keep devices away from children.
- Use only manufacturer-supplied and recommended attachments.
- Keep the device away from open flames and heat. Fire, electrical shock, damage, or injury may occur.
- Keep the pill dispenser on a flat, stable surface. Dropping may cause damage to the dispenser or pills. Flip it only if you're removing pills. The Internal pill container is removable for easy refilling.
- Water can damage the dispenser, batteries, or the medication inside. Do not immerse or expose this pill dispenser to water or other liquids. Keep in a cool, dry environment.
- Do not fill the internal pill container unless it is completely dry.
- The user, caregiver, or healthcare provider is solely responsible for organizing medications and assigning quantities and times.

- MOBI® is not liable for missed medications.
- If a dose is missed, contact your pharmacist or healthcare provider for instructions.
- Some medications should be kept in an airtight container. Ask your healthcare provider or pharmacist if you're unsure if this dispenser is right for you.
- Change the batteries immediately when a low battery symbol  appears in the upper right corner of the LCD. Never use mismatched, old, or damaged batteries. Never force the locking cover open.
- Do not mix Alkaline, Standard, or Rechargeable Batteries



## Battery Safety

- Use the batteries recommended.
- Make sure the battery compartment is secure.
- Fit batteries correctly, observing the proper polarities.
- Store unused batteries in their packaging and keep away from metal objects which could cause a short circuit.
- Remove dead batteries from the scales.
- Remove batteries from the scales if they will not be used for a prolonged period.

### Warnings:

- Never dispose of batteries in a fire.
- Never attempt to recharge ordinary batteries.
- Do not mix old and new batteries.
- Do not mix alkaline, standard (carbon-zinc), or rechargeable (ni-cad, ni-mg, etc.) batteries.

### Changing the battery:

You can access the battery by removing the door on scale's underside.

## Care & Maintenance

1. Wipe the exterior with a soft, damp cloth. Do not use harsh detergents.
2. If needed, remove the internal pill container, and gently rinse/wipe the individual pill compartments. Allow the container to dry fully before replacing it or refilling it.
3. Do not immerse the outer pill dispenser in water or any other liquid.
4. Make sure that the internal pill container is clean and dry before refilling it.
5. Please follow local and state guidelines for proper battery disposal.



## FCC Statement

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this device

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

FCC Radiation Exposure Statement:

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with a minimum distance of 20cm between the radiator & your body.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:



- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

## Limited Product Warranty

MOBI® warrants its products to be free from defects in material and workmanship for 90 days from the date of purchase. Should the product prove to be defective during the warranty period, we would either replace or repair the product at our option without charge. After the warranty period, a service charge will apply for replacing parts and labor costs. To obtain warranty service, you must return the product complete with its original parts & accessories and a copy of the original purchase, gift, or registry to us. This warranty does not cover damage caused by tampering, accident, misuse, abuse, improper maintenance, unauthorized modification, or connection to an improper power supply. There will be a charge for repair caused for such damage. This warranty excludes all incidental or consequential damages and any liability other than stated above. MOBI® reserves the right to modify or discontinue the offer by posting a notice on the app or website.



## Extended Manual

Scan to find the most up-to-date Instruction Manual & Warranty Details.



If you experience issues setting up your device, visit [www.mobiusa.com/gethelp](http://www.mobiusa.com/gethelp).

Email us at [service@mobiusa.com](mailto:service@mobiusa.com).

You can also reach us on **+1 (323) 904-3000**,  
Monday - Friday, 8:30 AM - 5:00 PM PST.



Visit [www.mobiusa.com](http://www.mobiusa.com)

©2024 MOBI® Technologies, Inc.  
All Rights Reserved.