



10 series Upper Arm Blood Pressure Monitor Model BP7450 Instruction Manual

Table of Contents

Introduction	2	8.3 Cleaning	43
Safety Instructions	2	8.4 Disposal	43
Intended Use	2	9. Optional Accessories	44
Receiving and Inspection	2	10. Specifications	45
Symbols Glossary	2	11. FCC Statement	49
Important Safety Information	3	12. Limited Warranty	50
1. Know Your Monitor	7	13. Guidance and Manufacturer's Declaration	52
1.1 Contents/Product Includes	7		
1.2 Monitor and Arm Cuff	7		
1.3 Display and Symbols	9		
1.4 Classification of BP (Blood Pressure)	13		
2. Preparation	14		
2.1 Installing Batteries	14		
2.2 Using the AC Adapter	15		
2.3 Setting the Date and Time	16		
2.4 Setting the TruRead Mode	18		
2.5 Blood Pressure Measurement Tips	20		
2.6 Applying the Arm Cuff	21		
2.7 Sitting Correctly	23		
2.8 Selecting a User ID	24		
3. Using Your Monitor without a Smart Device	25		
3.1 Taking a Measurement	25		
4. Using Your Monitor with a Smart Device	29		
4.1 Pairing Your Monitor with a Smart Device	29		
4.2 Taking a Measurement	30		
4.3 Disabling/Enabling Bluetooth	31		
5. Using Memory Function	32		
5.1 Viewing the Readings Stored in Memory	32		
5.2 Viewing the Average Reading	34		
6. Other Settings	35		
6.1 Deleting All Readings Stored in Memory	35		
6.2 Restoring Your Monitor to the Default Settings	36		
7. Error Messages and Troubleshooting	38		
8. Maintenance	42		
8.1 Maintenance	42		
8.2 Storage	42		

Introduction

Thank you for purchasing the OMRON BP7450 Blood Pressure Monitor. This monitor can be used by multiple users. This blood pressure monitor uses the oscillometric method of blood pressure measurement. This means this monitor detects your blood movement through your brachial artery and converts the movements into a digital reading.

Safety Instructions

This instruction manual provides you with important information about the OMRON BP7450 Blood Pressure Monitor. To ensure the safe and proper use of this monitor, READ and UNDERSTAND all of the safety and operating instructions. **If you do not understand these instructions or have any questions, contact 1-800-634-4350 before attempting to use this monitor. For specific information about your own blood pressure, consult with your physician.**

Intended Use

This device is a digital monitor intended for use in measuring blood pressure and pulse rate in adult patient population. The device detects the appearance of irregular heartbeats during measurement and gives a warning signal with readings.

Environments of Use: Home

Patient Population: Adult

Receiving and Inspection

Remove this monitor and other components from the packaging and inspect for damage. If this monitor or any other components is damaged, DO NOT USE and contact 1-800-634-4350.

Symbols Glossary

For symbol information, visit:

OmronHealthcare.com/symbols-glossary

Contraindications

- This monitor is contraindicated for use in ambulatory environments.
 - This monitor is contraindicated for use on aircraft.
-

Important Safety Information

Read the Important Safety Information in this instruction manual before using this monitor.

Follow this instruction manual thoroughly for your safety.

Keep for future reference. For specific information about your own blood pressure, CONSULT WITH YOUR PHYSICIAN.

Warning Indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury.

- DO NOT use this monitor on infants, toddlers, children or persons who cannot express themselves.
- DO NOT adjust medication based on readings from this blood pressure monitor. Take medication as prescribed by your physician. ONLY a physician is qualified to diagnose and treat high blood pressure.
- DO NOT use this monitor on an injured arm or an arm under medical treatment.
- DO NOT apply the arm cuff on your arm while on an intravenous drip or blood transfusion.
- DO NOT use this monitor in areas containing high frequency (HF) surgical equipment, magnetic resonance imaging (MRI) equipment, computerized tomography (CT) scanners. This may result in incorrect operation of the monitor and/or cause an inaccurate reading.
- DO NOT use this monitor in oxygen rich environments or near flammable gas.
- Consult with your physician before using this monitor if you have common arrhythmias such as atrial or ventricular premature beats or atrial fibrillation; arterial sclerosis; poor perfusion; diabetes; pregnancy; pre-eclampsia or renal disease. NOTE that any of these conditions in addition to patient motion, trembling, or shivering may affect the measurement reading.
- NEVER diagnose or treat yourself based on your readings. ALWAYS consult with your physician.
- To help avoid strangulation, keep the air tube and AC adapter cable away from infants, toddlers or children.
- This product contains small parts that may cause a choking hazard if swallowed by infants, toddlers or children.

Data Transmission

- This product emits radio frequencies (RF) in the 2.4 GHz band. DO NOT use this product in locations where RF is restricted, such as on an aircraft or in hospitals. Turn off the Bluetooth® feature in this monitor and remove batteries and/or unplug the AC adapter when in RF restricted areas. For further information on potential restrictions refer to documentation on the Bluetooth usage by the FCC.

AC Adapter Handling and Usage

- DO NOT use the AC adapter if this monitor or the AC adapter cable is damaged. If this monitor or the cable is damaged, turn off the power and unplug the AC adapter immediately.
- Plug the AC adapter into the appropriate voltage outlet. DO NOT use in a multi-outlet plug.
- NEVER plug in or unplug the AC adapter from the electric outlet with wet hands.
- DO NOT disassemble or attempt to repair the AC adapter.

Battery Handling and Usage

- Keep batteries out of the reach of infants, toddlers or children.
-



Caution

Indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury to the user or patient, or cause damage to the equipment or other property.

- Stop using this monitor and consult with your physician if you experience skin irritation or discomfort.
- Consult with your physician before using this monitor on an arm where intravascular access or therapy, or an arterio-venous (A-V) shunt, is present because of temporary interference to blood flow which could result in injury.
- Consult with your physician before using this monitor if you have had a mastectomy.
- Consult with your physician before using this monitor if you have severe blood flow problems or blood disorders as cuff inflation can cause bruising.
- DO NOT take measurements more often than necessary because bruising, due to blood flow interference, may occur.
- ONLY inflate the arm cuff when it is applied on your upper arm.
- Remove the arm cuff if it does not start deflating during a measurement.
- DO NOT use this monitor for any purpose other than measuring blood pressure.
- During measurement, make sure that no mobile device or any other electrical device that emits electromagnetic fields is within 12 inches (30 cm) of this monitor. This may result in incorrect operation of the monitor and/or cause an inaccurate reading.
- DO NOT disassemble or attempt to repair this monitor or other components. This may cause an inaccurate reading.
- DO NOT use in a location where there is moisture or a risk of water splashing this monitor. This may damage this monitor.
- DO NOT use this monitor in a moving vehicle such as in a car.
- DO NOT drop or subject this monitor to strong shocks or vibrations.

- **DO NOT** use this monitor in places with high or low humidity or high or low temperatures. Refer to section 10.
- During measurement, observe the arm to ensure that the monitor is not causing prolonged impairment to blood circulation.
- **DO NOT** use this monitor in high-use environments such as medical clinics or physician offices.
- **DO NOT** use this monitor with other medical electrical (ME) equipment simultaneously. This may result in incorrect operation of the monitor and/or cause an inaccurate reading.
- Avoid bathing, drinking alcohol or caffeine, smoking, exercising and eating for at least 30 minutes before taking a measurement.
- Rest for at least 5 minutes before taking a measurement.
- Remove tight-fitting, thick clothing from your arm while taking a measurement.
- Remain still and **DO NOT** talk while taking a measurement.
- **ONLY** use the arm cuff on persons whose arm circumference is within the specified range of the cuff.
- Ensure that this monitor has acclimated to room temperature before taking a measurement. Taking a measurement after an extreme temperature change could lead to an inaccurate reading. OMRON recommends waiting for approximately 2 hours for the monitor to warm up or cool down when the monitor is used in an environment within the temperature specified as operating conditions after it is stored either at the maximum or at the minimum storage temperature. For additional information of operating and storage/transport temperature, refer to section 10.
- **DO NOT** use this monitor after the durable period has ended. Refer to section 10.
- **DO NOT** crease the arm cuff or the air tube excessively.
- **DO NOT** fold or kink the air tube while taking a measurement. This may cause an injury by interrupting blood flow.
- To unplug the air plug, pull on the plastic air plug at the base of the tube, not the tube itself.
- **ONLY** use the AC adapter, arm cuff, batteries and accessories specified for this monitor. Use of unsupported AC adapters, arm cuffs and batteries may damage and/or may be hazardous to this monitor.
- **ONLY** use the approved arm cuff for this monitor. Use of other arm cuffs may result in incorrect readings.
- Inflating to a higher pressure than necessary may result in bruising of the arm where the cuff is applied. **NOTE:** refer to “If your systolic pressure is more than 210 mmHg” in section 3 of this instruction manual for additional information.

Data Transmission

- DO NOT replace batteries or unplug the AC adapter while your readings are being transferred to your smart device. This may result in incorrect operation of this monitor and failure to transfer your blood pressure data.

AC Adapter Handling and Usage

- Fully insert the AC adapter into the outlet.
- When unplugging the AC adapter from the outlet, be sure to safely pull from the AC adapter. DO NOT pull from the AC adapter cable.
- When handling the AC adapter cable:
 - DO NOT damage it.
 - DO NOT break it.
 - DO NOT tamper with it.
 - DO NOT forcibly bend or pull it.
 - DO NOT twist it.
 - DO NOT use it if it is gathered in a bundle.
 - DO NOT pinch it.
 - DO NOT place it under heavy objects.
- Wipe any dust off of the AC adapter.
- Unplug the AC adapter when not in use.
- Unplug the AC adapter before cleaning this monitor.

Battery Handling and Usage

- DO NOT insert batteries with their polarities incorrectly aligned.
 - ONLY use 4 “AA” alkaline or manganese batteries with this monitor. DO NOT use other types of batteries. DO NOT use new and used batteries together. DO NOT use different brands of batteries together.
 - Remove batteries if this monitor will not be used for a long period of time.
 - If battery fluid should get in your eyes, immediately rinse with plenty of clean water. Consult with your physician immediately.
 - If battery fluid should get on your skin, wash your skin immediately with plenty of clean, lukewarm water. If irritation, injury or pain persists, consult with your physician.
 - DO NOT use batteries after their expiration date.
 - Periodically check batteries to ensure they are in good working condition.
-


1. Know Your Monitor

1.1 Contents/Product Includes

Monitor, Easy Wrap ComFit Cuff, AC adapter (HEM-ADPTW5), 4 “AA” batteries, instruction manual, quick start guide


1.2 Monitor and Arm Cuff

Front of the monitor

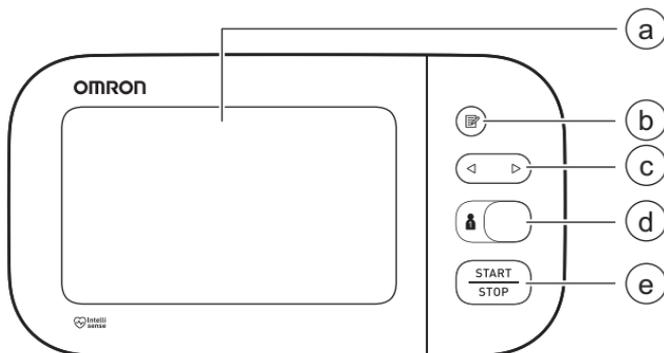

- (a) Display
- (b) [Memory] button
- (c) [Forward/Backward] buttons
- (d) [User ID selection] switch
- (e) [START/STOP] button


Left and right side of the monitor

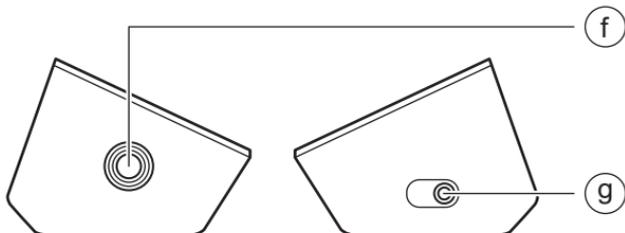

- (f) Air jack
- (g) AC adapter jack

1. Know Your Monitor

Back and bottom of the monitor

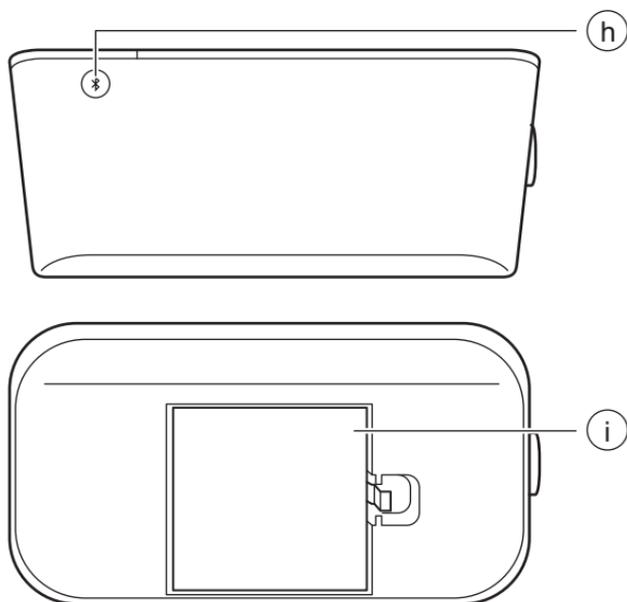

Ⓜ [Connection] button
(Date/Time)


Ⓜ Battery compartment

Arm cuff


(Upper arm circumference: 9" - 17" (22 - 42 cm))

Ⓜ Air plug

Ⓜ Air tube

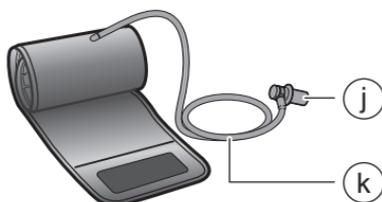

1.3 Display and Symbols

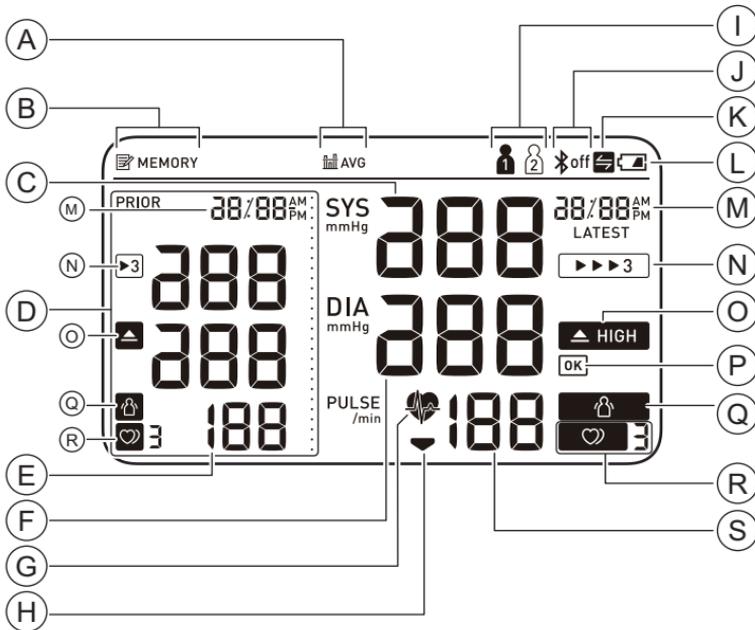

- | | |
|--------------------------------------|--|
| (A) Average value symbol | (K) Sync symbol |
| (B) Memory symbol | (L) Low battery symbol / Depleted battery symbol |
| (C) Systolic blood pressure reading | (M) Date/time display |
| (D) Comparison display | (N) TruRead symbol |
| (E) Pulse display / Memory number | (O) High symbol |
| (F) Diastolic blood pressure reading | (P) OK symbol |
| (G) Heartbeat symbol | (Q) Movement error symbol |
| (H) Deflation symbol | (R) Irregular heartbeat symbol |
| (I) User ID symbol | (S) Pulse display |
| (J) Bluetooth ON / OFF symbol | |

1. Know Your Monitor

(A)	 AVG	Average value symbol Flashes when viewing the average of the latest 2 or 3 readings taken within the most recent 10 minute timeframe.
(B)	 MEMORY	Memory symbol Appears when viewing readings stored in the memory.
(C)	Systolic blood pressure reading	
(D)	Comparison display Displays past readings. Refer to section 5 for more information.	
(E)	Pulse display / Memory number When pressing the  button, the memory number appears for approximately 1.5 seconds before the pulse rate appears on the display.	
(F)	Diastolic blood pressure reading	
(G)		Heartbeat symbol Flashes while taking a measurement.
(H)		Deflation symbol Appears during cuff deflation.
(I)		User ID symbol The user ID number appears when selecting the user with [user ID selection] switch.
(J)		Bluetooth ON symbol Appears when your readings are being transferred.
	 off	Bluetooth OFF symbol Appears when Bluetooth is disabled.

(K)		<p>Sync symbol Flashes/appears when your data needs to be transferred because the stored memory is either almost, or completely full. Once you pair your monitor with your smart device, transfer your readings immediately before the monitor deletes the oldest reading. Up to 100 readings can be stored in the internal memory of your monitor per user.</p>
(L)		<p>Low battery symbol Flashes when batteries are low.</p>
		<p>Depleted battery symbol Appears when batteries are depleted.</p>
(M)	Date/time display	
(N)		<p>TruRead symbol Appears when a measurement is taken if the TruRead mode is enabled. The TruRead mode automatically takes 3 consecutive readings at intervals you can customize and displays the average. To enable the TruRead mode, refer to sub-section 2.4. To use the TruRead mode, refer to sub-section 3.1. To view your readings in TruRead measurement, refer to sub-section 5.1.</p>
(O)		<p>High symbol Appears if your systolic blood pressure is ≥ 130 mmHg or the diastolic blood pressure is ≥ 80 mmHg. Refer to sub-section 1.4 for more information.</p>
(P)		<p>OK symbol Flashes when your monitor is connected to your smart device or readings are transferred successfully.</p>

1. Know Your Monitor

Q		<p>Movement error symbol Appears with your reading when your body is moving during a measurement. If it appears, remove the arm cuff and wait for 2 to 3 minutes. Then, reapply the cuff, remain still and take another measurement.</p>
R	<p></p> <p>or</p> <p></p>	<p>Irregular heartbeat symbol Appears when irregular rhythm is detected 2 or more times during a measurement. The symbol with the number appears as below:</p> <p>Irregular heartbeat is detected 2 times:  2 Irregular heartbeat is detected 3 times or more:  3</p> <p>An irregular heartbeat rhythm is defined as a rhythm that is 25% less or 25% more than the average rhythm detected while your monitor is measuring blood pressure.</p> <p></p> <p>Irregular heartbeat symbol may be displayed when the monitor detects arm movements. Refer to section 7 for more information.</p> <p>If it continues to appear, we recommend you to consult with and follow the directions of your physician.</p>
S		<p>Pulse display Pulse rate appears after the measurement.</p>

1.4 Classification of BP (Blood Pressure)

ACC/AHA 2017 Hypertension guidelines define categories of BP in adults as below.

Categories of BP in Adults*

BP Category	Systolic BP	Diastolic BP
Normal	< 120 mmHg	and < 80 mmHg
Elevated	120 - 129 mmHg	and < 80 mmHg
Hypertension		
Stage 1	130 - 139 mmHg	or 80 - 89 mmHg
Stage 2	≥ 140 mmHg	or ≥ 90 mmHg

* Individuals with Systolic BP and Diastolic BP in 2 categories should be designated to the higher BP category.

BP indicates blood pressure (based on an average of ≥ 2 careful readings obtained on ≥ 2 occasions).

Source: ACC/AHA 2017 High Blood Pressure Clinical Practice Guideline


⚠Warning

- NEVER diagnose or treat yourself based on your readings. ALWAYS consult with your physician.

2. Preparation

2.1 Installing Batteries

1. Push down the hook of the battery cover and pull downward.

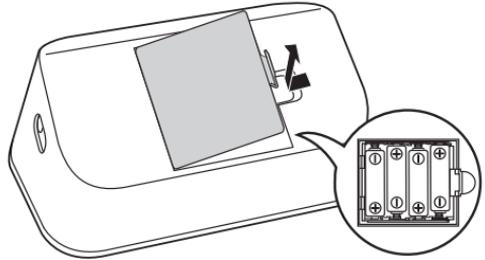

2. Insert 4 “AA” batteries as indicated in the battery compartment.
3. Close the battery cover.

As soon as batteries are correctly installed, the year flashes on the display.

Refer to sub-section 2.3 to set the date and time manually, or refer to sub-section 4.1 to pair your monitor to set the date and time automatically.



Note


- When the “” symbol flashes on the display, it is recommended you replace batteries.
- To replace batteries, turn your monitor off and remove all batteries. Then replace with 4 new batteries at the same time.
- Replacing batteries will not delete previous readings.
- The supplied batteries may have a shorter life span than new batteries.
- Disposal of used batteries should be carried out in accordance with local regulations.

2.2 Using the AC Adapter

Note

- Make sure not to place your monitor in a location where it is difficult to plug and unplug the AC adapter.
- We recommend keeping batteries in your monitor at all times, even if you choose to use the AC adapter. If only the AC adapter is used without keeping batteries in your monitor, you may need to reset the date and time each time you unplug and plug back the AC adapter. The readings will not be deleted.

1. Insert the AC adapter plug into the AC adapter jack on the right side of your monitor.

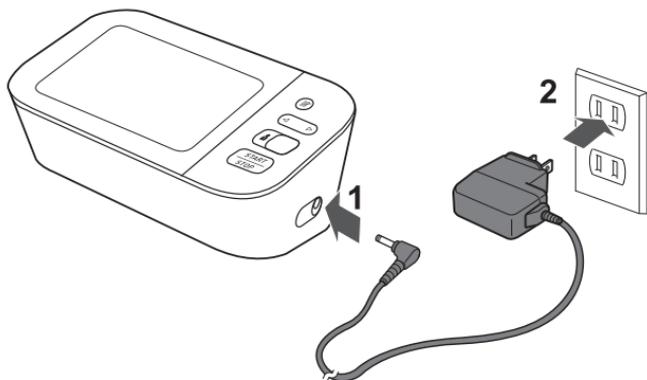

2. Plug the AC adapter into an electrical outlet.

To unplug the AC adapter, unplug the AC adapter from the electrical outlet and then remove the AC adapter plug from the monitor.

Caution

- Only use the AC adapter, the arm cuff, batteries and accessories specified for this monitor. Use of unsupported AC adapters, arm cuffs and batteries may damage and/or may be hazardous to this monitor.
-


2.3 Setting the Date and Time

Note

- Set your monitor to the correct date and time before taking a measurement for the first time.
- The date and time will automatically be set when the monitor is paired with “OMRON connect US/CAN” app. Refer to sub-section 4.1.
- As soon as batteries are correctly installed, the year flashes on the display of your monitor.

1. If nothing is displayed, press  button (rear surface of the monitor) to display the year.

If the year is already displayed, skip this step.

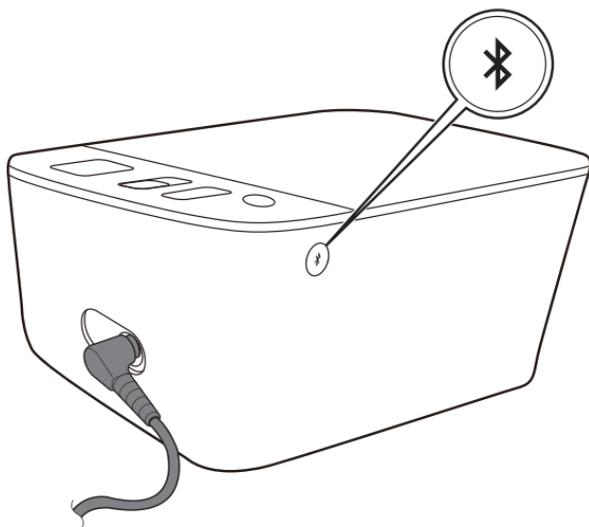


2. Press the  or  button to adjust the year.

Press and hold the  button to quickly advance the years.

Press and hold the  button to quickly go backwards through the years.

-
3. Press the  button to confirm the year, then the month will flash. Repeat the same steps to adjust the month, day, hour, and minute.

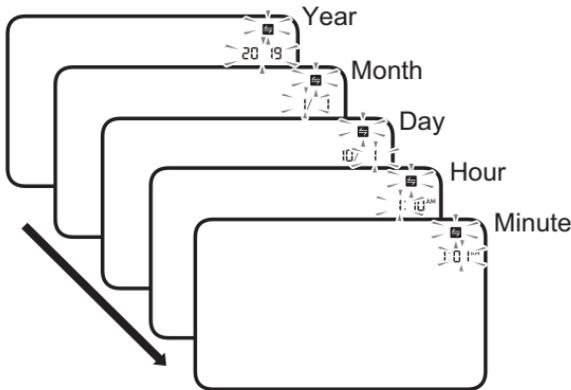

-
4. Press the  button to store the setting.

TruRead mode setting is displayed 3 seconds after the minutes have been set in this sequence. Refer to sub-section 2.4.

Note

- If you need to reset the date and time, replace batteries or press the  button, then set the date and time as shown in above steps.
- If the date and time are not set, “-:--” flashes during or after measurement.

2.4 Setting the TruRead Mode

The TruRead mode takes 3 consecutive measurements. The monitor will inflate, take a measurement, and deflate - 3 times, separated by a short interval between each measurement.

The TruRead mode is set to “oFF” by default.

1. Set or bypass the date and time settings.

Follow the steps in sub-section 2.3 and set or bypass the date and time settings.


After it is completed, the “▶▶▶3” symbol appears on the display.



2. To enable TruRead mode, press the ◀ or ▶ button.

“oN” will appear in the display.

To keep TruRead mode disabled, press the Ⓢ button (rear surface of the monitor).

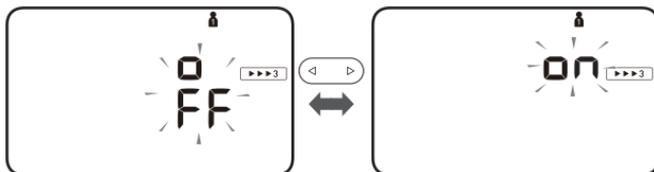


3. Press the Ⓢ button to confirm.

If “oN” is selected, the interval setting appears.



4. Press the ◀ or ▶ button to change the interval.

It can be set 15, 30, 60, or 120 seconds.

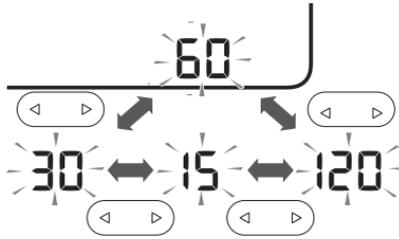

5. Press the Ⓢ button to confirm.

Your monitor will automatically turn off in 3 seconds.



To take a measurement in TruRead mode, refer to “Using the TruRead Mode” in sub-section 3.1.

2.5 Blood Pressure Measurement Tips

To help ensure an accurate reading, follow these directions:

- Stress raises blood pressure. Avoid taking measurements during stressful times.
- Measurements should be taken in a quiet place.
- It is important to take measurements at the same times each day. Taking measurements in the morning and in the evening is recommended.
- Remember to have a record of your blood pressure and pulse readings for your physician. A single measurement does not provide an accurate indication of your true blood pressure. You need to take and record several readings over a period of time.

Caution

- **Avoid bathing, drinking alcohol or caffeine, smoking, exercising and eating for at least 30 minutes before taking a measurement.**
 - **Rest for at least 5 minutes before taking a measurement.**
-

2.6 Applying the Arm Cuff


Note

- The following steps are for applying the arm cuff to your left arm. When you take a measurement on your right arm, follow the cuff wrap instructions for a right arm use mentioned under “To take a measurement on your right arm” located at the end of this sub-section.
- The blood pressure can differ between the right arm and the left arm, and the measured blood pressure readings can be different. OMRON recommends to always use the same arm for measurement. If the readings between both arms differ substantially, check with your physician to determine as to which arm to use for your measurements.


⚠ Caution

- Remove tight-fitting, thick clothing from your arm while taking a measurement.

1. **Plug the arm cuff into your monitor by inserting the air plug into the air jack securely until it clicks.**

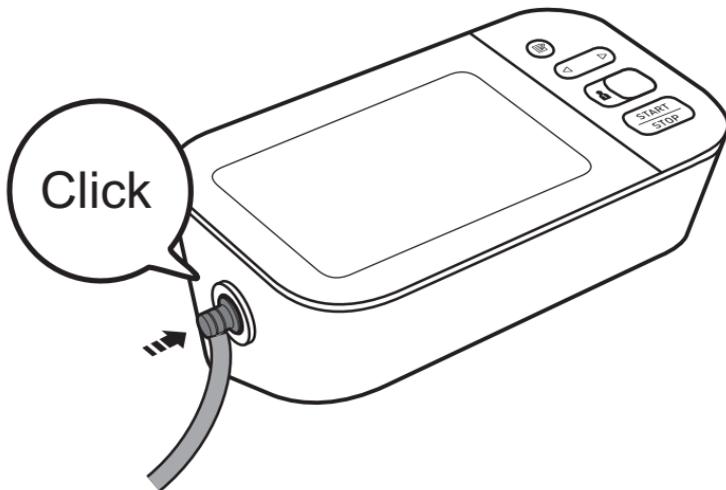


2. Apply the arm cuff to your left upper arm.

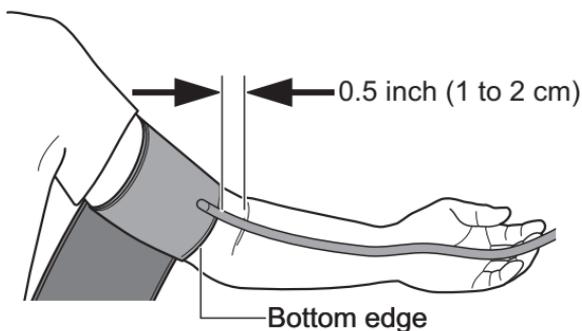

Note

- The bottom edge of the arm cuff should be 0.5 inch (1 to 2 cm) above the inside elbow. The air tube is on the inside of your arm and aligned with your middle finger.

3. Make sure that the air tube is positioned on the inside of your arm and wrap the cuff securely, so it can not move around your arm.

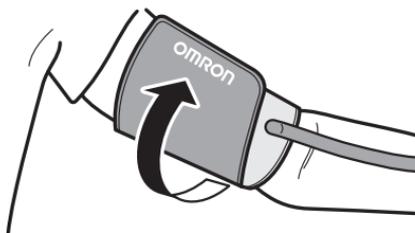

To take a measurement on your right arm

When you take a measurement on your right arm, the air tube should run along the side of your elbow, along the bottom of your arm. Be careful not to rest your arm on the air tube.



2.7 Sitting Correctly

To take a measurement, you need to be relaxed and comfortably seated in a room at a comfortable temperature. Place your arm on the table.


- Sit in a comfortable chair with your back and arm supported.
- Keep your feet flat and your legs uncrossed.
- The arm cuff should be placed on your arm at the same level as your heart, with the arm resting comfortably on a table.



2.8 Selecting a User ID

Your monitor is designed to take measurements and store the readings in the memory for 2 people using the [user ID selection] switch.

1. Select your user ID (1 or 2).

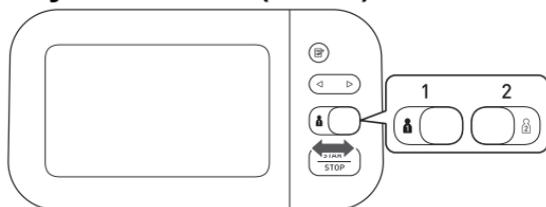

3. Using Your Monitor without a Smart Device

3.1 Taking a Measurement

You can use your monitor without pairing to a smart device. To pair your monitor with a smart device, refer to sub-section 4.1.

Note

- To stop the measurement, press the [START/STOP] button once to deflate the arm cuff.

⚠ Caution

- **DO NOT use this monitor with other medical electrical (ME) equipment simultaneously. This may result in incorrect operation of the monitor and/or cause an inaccurate reading.**
- **Remain still and DO NOT talk while taking a measurement.**

1. Press the [START/STOP] button.


All symbols appear on the display before starting the measurement.

2. Remain still and do not talk until the entire measurement process is complete.

The “” symbol flashes at every heartbeat.

After your monitor has detected your blood pressure and pulse rate, the cuff automatically deflates. Your blood pressure and pulse rate are displayed with the previous reading.

If either the systolic or the diastolic reading is high (refer to sub-section 1.3), the “ HIGH” symbol appears.

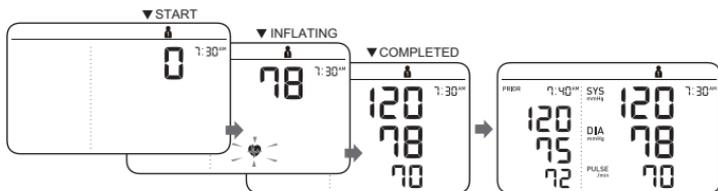

3. Press the [START/STOP] button to turn your monitor off.

Note

- Your monitor automatically turns off after 3 minutes.
- Wait 2 - 3 minutes between measurements. The wait time allows the arteries to decompress and return to their pre-measurement state. You may need to increase the wait time depending on your individual physiological characteristics.

■ Using the TruRead Mode

Set the TruRead mode “on” to use the TruRead mode. Refer to subsection 2.4.


1. Select your user ID.

2. Press the [START/STOP] button.

The arm cuff starts to inflate automatically.

After the first measurement is complete, wait for the next measurement to start.


The interval you selected during TruRead settings will appear on the display.

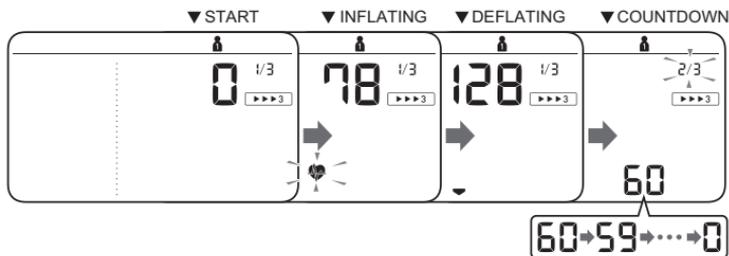


3. Using Your Monitor without a Smart Device

After the second measurement is complete, wait for the next measurement to start.

The interval you selected during TruRead settings will appear on the display.

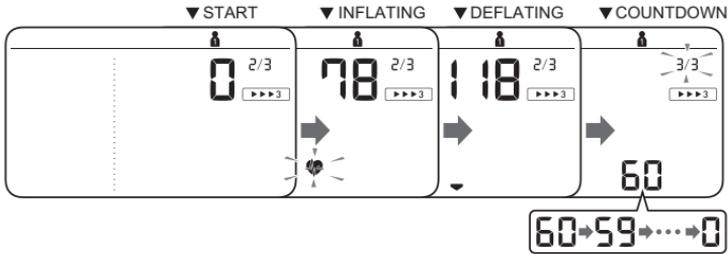

After the third measurement is complete, the average for the 3 measurements appears on the display.

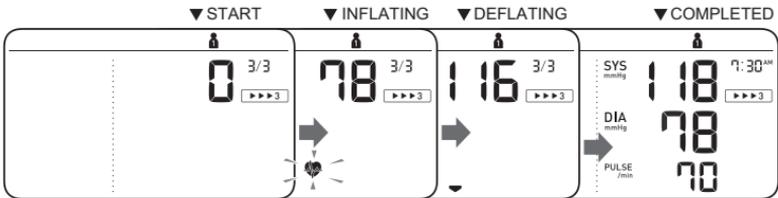

To view readings in TruRead measurement, refer to subsection 5.1.

3. Using Your Monitor without a Smart Device

■ Using the Guest Mode

Your monitor stores readings for 2 users in the memory. The guest mode can be used to take a single measurement for another user.

Note

- No readings are stored in the memory, and the TruRead mode is not available when the guest mode is selected.

1. **Press and hold the [START/STOP] button for more than 3 seconds.**



The user ID symbol and the date/time display will disappear.

2. **Release the [START/STOP] button when the user ID and date/time display turn off.**

The arm cuff automatically starts to inflate.

■ If your systolic pressure is more than 210 mmHg

If your systolic pressure is more than 210 mmHg, take a measurement manually.

After the arm cuff starts to inflate, press and hold the [START/STOP] button until the monitor inflates 30 to 40 mmHg higher than your expected systolic pressure.

Note

- Do not inflate above 299 mmHg.

⚠ Caution

- Inflating to a higher pressure than necessary may result in bruising where the arm cuff is applied.

4. Using Your Monitor with a Smart Device

4.1 Pairing Your Monitor with a Smart Device

The date and time on your monitor will automatically be set when you pair your monitor with your smart device.

Review the list of compatible smart devices at OmronHealthcare.com


1. **Enable Bluetooth on your smart device.**
2. **Download and install the free “OMRON connect US/CAN” app onto your smart device.**



If you already have the “OMRON connect US/CAN” app and have created your account, open the app and add your new monitor.

3. **Open the app and follow the pairing instructions shown on your smart device.**
4. **Confirm that your monitor is connected successfully.**


When your monitor is connected successfully to your smart device, the “OK” symbol appears.

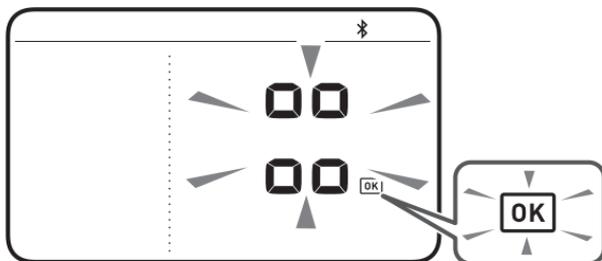

5. Press [START/STOP] button to turn your monitor off.

Note

- If the “Err” symbol appears, follow the instructions in the “OMRON connect US/CAN” app.

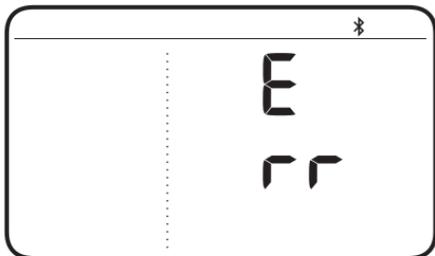

- Your monitor automatically turns off after 2 minutes.
- Be aware that OMRON will not be responsible for the loss of data and/or information in the app.
- “OMRON connect US/CAN” is the only app that we recommend to use with your monitor to transfer data correctly.

4.2 Taking a Measurement

Refer to sub-section 3.1 to start a measurement. Open the app on your smart device and follow the instructions.

Note

- Make sure that your monitor is within 16 ft. (5m) of your smart device.

■ Transferring Data Automatically/Manually

- The readings will be transferred to your smart device within 1 hour after taking a measurement.
- To transfer the readings manually, press the  button.

4.3 Disabling/Enabling Bluetooth


⚠Warning

- This product emits radio frequencies (RF) in the 2.4 GHz band. **DO NOT** use this product in locations where RF is restricted, such as on an aircraft or in hospitals. Turn off the Bluetooth feature in this monitor and remove batteries when in RF restricted areas. For further information on potential restrictions refer to documentation on the Bluetooth usage by the FCC.

Bluetooth is enabled as default.

1. **When your monitor is off, press and hold the  button (rear surface of the monitor) for more than 10 seconds.**


“oFF” flashes 3 times.

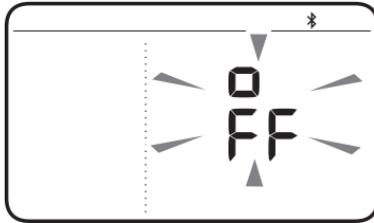

2. **Press the [START/STOP] button to turn your monitor off.**

Note

- To enable Bluetooth, press and hold the  button for more than 2 seconds. “on” appears on the display.
- Your monitor automatically turns off 2 minutes after the completion.
- When Bluetooth is disabled, the “ off” symbol appears.

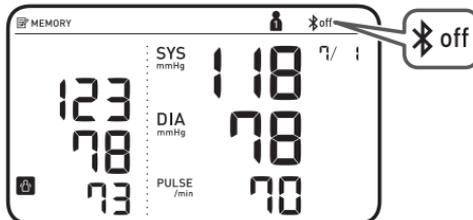

5. Using Memory Function

Your monitor automatically stores up to 100 readings per user.


Note

- This monitor can be used by multiple users. Refer to sub-section 2.8 to switch user. Make sure which user ID stores your own readings.

5.1 Viewing the Readings Stored in Memory

1. Select your user ID.

2. Press the button.

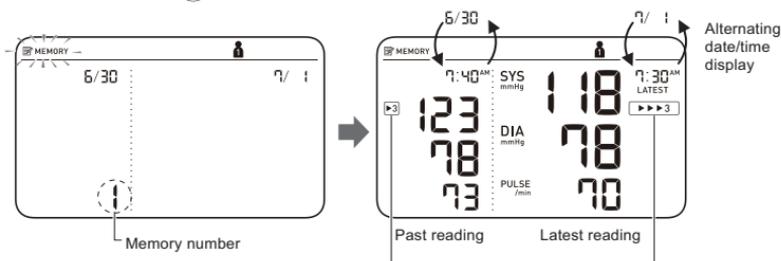


If the past reading was taken in TruRead mode, the “” symbol appears and the reading shown on the left is the average of 3 measurements.

If the latest reading was taken in TruRead mode, the “” symbol appears and the reading shown on the right is the average of 3 measurements.

After a memory number appears, the latest reading will appear with a past reading. A past reading will appear in the comparison display on the left.

Note

- If the reading is high (refer to sub-section 1.3), the “ HIGH” symbol appears.
- If there are no readings stored in the memory, the screen to the right is displayed.
- If the memory is full, new readings will replace the old readings, starting with the oldest reading.

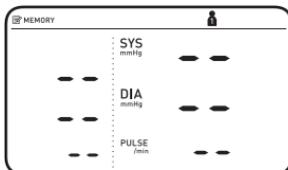

3. Press the ◀ or ▶ button repeatedly to scroll through the previous readings stored in the memory.

◀ : To view the older readings

▶ : To view the more recent readings


To view individual readings in TruRead measurement

The readings with the TruRead symbols (▶3 / ▶▶▶3) are the averages for the consecutive 3 measurements. To view the individual readings, press the Ⓢ button (rear surface of the monitor) while the TruRead average is displayed.

1. When a TruRead reading is the past reading displayed on the left side:



2. When a TruRead reading is the latest reading displayed on the right side:

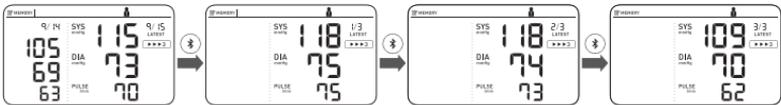

3. When TruRead readings are the past and the latest readings displayed in both of left and right side:



5.2 Viewing the Average Reading


Your monitor can calculate an average reading based on the latest 2 or 3 readings taken within the most recent 10 minute timeframe.

1. **Select your user ID.**
2. **When your monitor is off, press and hold the  button for more than 3 seconds.**

The “ AVG” symbol flashes.

Note

- If there are only 2 readings in the memory for that period, the average will be based on these 2 readings.
- If the average reading is high (refer to sub-section 1.3), the “” symbol appears.

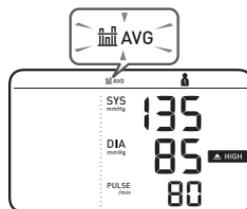

6. Other Settings

6.1 Deleting All Readings Stored in Memory


1. Select your user ID.

2. Press the button.

After a memory number appears, the latest reading will appear with a past reading.


3. While holding the button down, press and hold the [START/STOP] button for more than 3 seconds.

The display to the right flashes.

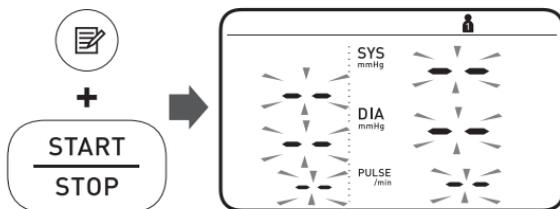

4. While holding the button down, press and hold the [START/STOP] button for more than 3 seconds.

The display stops flashing, then all readings are deleted.

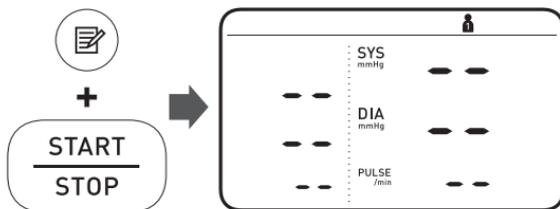

5. Press [START/STOP] button to turn your monitor off.

Note


- All readings will be deleted. You cannot partially delete the readings stored in the memory.
- Your monitor automatically turns off after 2 minutes.

6.2 Restoring Your Monitor to the Default Settings

To delete all the information stored in your monitor, follow the instructions below. Make sure that your monitor is turned off.


-
1. While holding the  button (rear surface of the monitor) down, press and hold the [START/STOP] button for more than 3 seconds.

The display to the right flashes.

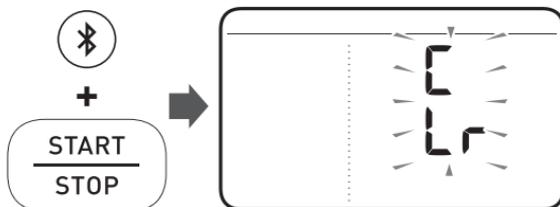

-
2. While holding the  button down, press and hold the [START/STOP] button for more than 3 seconds.

The display stops flashing, then your monitor is restored to the default settings.

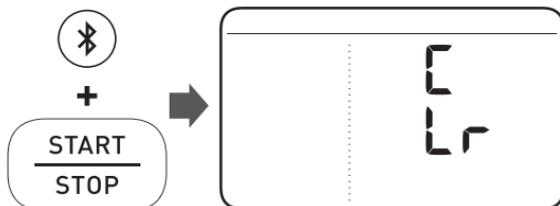

3. Press [START/STOP] button to turn your monitor off.

Note

- Reverting to the default setting of your monitor does not delete the information in the app.
- Your monitor automatically turns off after 3 minutes.
- You will need to re-pair your monitor when using the monitor again. Without re-pairing, your readings will not be transferred to your app.

7. Error Messages and Troubleshooting

If any of the below problems occur during measurement, check to make sure that no other electrical device is within 12 inches (30 cm). If the problem persists, refer to the table below.

Display/ Problem	Possible Cause	Solution
E1 appears or the arm cuff does not inflate.	Air plug is not completely plugged into the monitor.	Insert the air plug securely. Refer to sub-section 2.6.
	The arm cuff is not applied correctly.	Apply the arm cuff correctly, then take another measurement. Refer to sub-section 2.6.
	Air is leaking from the arm cuff.	Replace the arm cuff with a new one. Refer to section 9.
E2 appears or a measurement cannot be complete after the arm cuff inflates.	You move or talk during a measurement and the arm cuff does not inflate sufficiently.	Remain still and do not talk during a measurement. If "E2" appears repeatedly, inflate the arm cuff manually until the systolic pressure is 30 to 40 mmHg above your previous readings. Refer to sub-section 3.1.
	The systolic pressure is above 210 mmHg and a measurement cannot be taken.	Inflate the arm cuff manually until the systolic pressure is 30 to 40 mmHg above your previous readings. Refer to sub-section 3.1.

Display/ Problem	Possible Cause	Solution
 appears	The arm cuff is inflated exceeding the maximum allowable pressure.	Do not touch the arm cuff and/or bend the air tube while taking a measurement. If inflating the arm cuff manually, refer to the end of sub-section 3.1.
 appears	You move or talk during a measurement. Vibrations disrupt a measurement.	Remain still and do not talk during a measurement.
 appears  /  3 appears  does not flash during a measurement	The pulse rate is not detected correctly.	Apply the arm cuff correctly, then take another measurement. Refer to sub-section 2.6. Remain still and sit correctly during a measurement. If the “  3” symbol continues to appear, we recommend you to consult with your physician.
 appears	The monitor has malfunctioned.	Press the [START/STOP] button again. If “Er” still appears, contact customer service. Refer to section 12.

7. Error Messages and Troubleshooting

Display/ Problem	Possible Cause	Solution
 appears	The monitor cannot connect to a smart device or transmit data correctly.	Follow the instructions shown in the “OMRON connect US/CAN” app. If the “Err” symbol still appears after checking the app, contact customer service. Refer to section 12.
 flashes	The monitor is waiting for pairing with the smart device.	Refer to sub-section 4.1 for pairing your monitor with your smart device, or press [START/STOP] button to cancel pairing and turn your monitor off.
 flashes	The monitor is ready to transfer your readings to the smart device.	Open the “OMRON connect US/CAN” app to transfer your readings.
 flashes	There are more than 80 readings in memory to be transferred.	Pair or transfer your readings to the “OMRON connect US/CAN” app so you can keep them in memory in the app, and this error symbol disappears.
	The date and time is not set.	
 appears	There are 100 readings in memory to be transferred.	
 flashes	Batteries are low.	Replacing all 4 batteries with new ones is recommended. Refer to sub-section 2.1.

Display/ Problem	Possible Cause	Solution
 appears or the monitor is turned off unexpectedly during a measurement.	Batteries are depleted.	Immediately replace all 4 batteries with new ones. Refer to sub-section 2.1.
Nothing appears on the display of the monitor.	Batteries are depleted.	Check the battery installation for proper placement. Refer to sub-section 2.1.
Readings appear too high or too low.	Blood pressure varies constantly. Many factors including stress, time of day, and/or how you apply the arm cuff, may affect your blood pressure. Review sub-sections 2.5 - 2.7 and section 3.	
Any other communication issue occurs.	Follow the instructions shown in the smart device, or visit the "Support" section in the "OMRON connect US/CAN" app for further help. If the problem still persists, contact customer service. Refer to section 12.	
Any other problems occur.	Press the [START/STOP] button to turn the monitor off, then press it again to take a measurement. If the problem continues, remove all batteries and wait for 30 seconds. Then re-install batteries. If the problem still persists, contact customer service. Refer to section 12.	

8. Maintenance

8.1 Maintenance

To protect your monitor from damage, follow the directions below:

- Changes or modifications not approved by the manufacturer will void the user warranty.

⚠ Caution

- **DO NOT** disassemble or attempt to repair this monitor or other components. This may cause an inaccurate reading.

8.2 Storage


Store your monitor and other components in a clean, safe location.

1. Remove the arm cuff from the monitor.

⚠ Caution

- To unplug the air plug, pull on the plastic air plug at the base of the tube, not the tube itself.

2. Gently fold the air tube into the arm cuff.

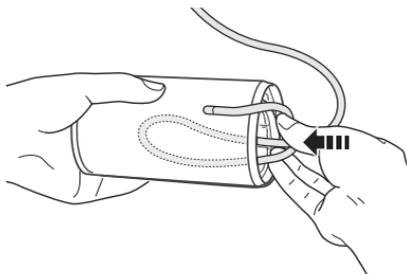

Note

- Do not bend or crease the air tube excessively.

Do not store your monitor and other components:

- If your monitor and other components are wet.
- In locations exposed to extreme temperatures, humidity, direct sunlight, dust or corrosive vapors such as bleach.
- In locations exposed to vibrations or shocks.

8.3 Cleaning

- Do not use any abrasive or volatile cleaners.
- Use a soft dry cloth or a soft cloth moistened with mild (neutral) detergent to clean your monitor and arm cuff and then wipe them with a dry cloth.
- Do not wash or immerse your monitor and arm cuff or other components in water.
- Do not use gasoline, thinners or similar solvents to clean your monitor and arm cuff or other components.

8.4 Disposal


Dispose of your monitor, other components and optional accessories according to applicable local regulations.

Unlawful disposal may cause environmental pollution.

9. Optional Accessories


Easy Wrap ComFit Cuff

Arm circumference
9" - 17" (22 - 42 cm)

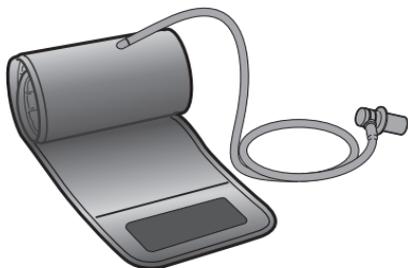


Model: HEM-FL31-B REF HEM-FL31-BZ

AC Adapter

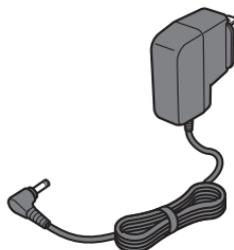

Model: HEM-ADPTW5

Protective LCD Cover

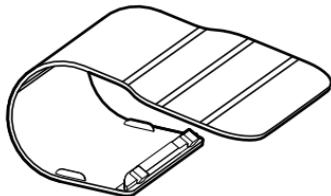

Model: HEM-CACO-734

For further information, visit our
website at OmronHealthcare.com

Note

- Do not throw the air plug away. The air plug can be applicable to the optional cuff.

10. Specifications

Model	BP7450 <small>REF</small> HEM-7342T-Z
Display	LCD digital display
Cuff pressure range	0 to 299 mmHg
Blood pressure measurement range	SYS: 60 to 260 mmHg DIA: 40 to 215 mmHg
Pulse measurement range	40 to 180 beats / min.
Accuracy	Pressure: ± 3 mmHg Pulse: $\pm 5\%$ of display reading
Inflation	Automatic by electric pump
Deflation	Automatic pressure release valve
Measurement method	Oscillometric method
Transmission method	Bluetooth [®] low energy technology
Wireless communication	Frequency range: 2.4 GHz (2400 - 2483.5 MHz) Modulation: GFSK Effective radiated power: < 20 dBm
Operation mode	Continuous operation
IP classification	Monitor: IP 20 AC adapter: IP21
Rating	DC6 V 4.0 W
Power source	4 "AA" batteries 1.5 V or AC adapter (INPUT AC100 - 240 V 50 - 60 Hz 0.12 - 0.065 A)
Battery life	Approximately 700 measurements (using new alkaline batteries)
Durable period (Service life)	Monitor: 5 years Arm cuff: 5 years AC adapter: 5 years
Operating conditions	+50°F to +104°F (+10°C to +40°C) / 15 to 90% RH (non-condensing) / 800 to 1060 hPa
Storage / Transport conditions	-4°F to +140°F (-20°C to +60°C) / 10 to 90% RH (non-condensing)

10. Specifications

Weight	Monitor: Approximately 16.6 oz. (470 g) not including batteries Arm cuff: Approximately 5.8 oz. (163 g)
Dimensions	Monitor: Approximately 7.5" (w) × 3.1" (h) × 4.1" (l) (191 mm × 79 mm × 105 mm) Arm cuff: Approximately 5.7" × 20.9" (air tube: 29.5") (145 mm × 532 mm (air tube: 750 mm))
Cuff circumference applicable to the monitor	9" to 17" (22 to 42 cm)
Memory	Stores up to 100 readings per user
Contents	Monitor, Easy Wrap ComFit Cuff, AC adapter (HEM-ADPTW5), 4 "AA" batteries, instruction manual, quick start guide
Protection against electric shock	Internally powered ME equipment (When using only batteries) Class II ME equipment (AC adapter)
Applied part	Type BF (arm cuff)

Note

- These specifications are subject to change without notice.
- In the clinical validation study, K5 was used on 85 subjects for determination of diastolic blood pressure.
- This monitor complies with the requirements of ISO 81060-2:2013.
- This monitor has not been validated for use on pregnant patients.
- IP classification is degrees of protection provided by enclosures in accordance with IEC 60529. This monitor and AC adapter are protected against solid foreign objects of 0.5 inches (12.5 mm) diameter and greater such as a finger. The AC adapter is protected against vertically falling water drops which may cause issues during a normal operation.

About a wireless communication interference

This product operates in an unlicensed ISM band at 2.4 GHz. In the event this product is used near other wireless devices such as microwave and wireless LAN, which operate on the same frequency band as this product, there is a possibility that interference may occur.

If interference occurs, stop the operation of the other devices or relocate this product away from other wireless devices before attempting to use it.

Pictogram on the Product

 ART.	Types of marks used to indicate the correct positioning of the artery for the upper arm
	Identifies the type of cuff that is compatible with the blood pressure monitor

11. FCC Statement

FCC CAUTION

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note

- This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
 - Reorient or relocate the receiving antenna.
 - Increase the separation between the equipment and receiver.
 - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 - Consult the dealer or an experienced radio/TV technician for help.

This transmitter must not be co-located or operated in conjunction with any other antenna or transmitter.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment and meets the FCC radio frequency (RF) Exposure Guidelines. This equipment has very low levels of RF energy that are deemed to comply without testing of specific absorption ratio (SAR).

12. Limited Warranty

Your BP7450 Blood Pressure Monitor, excluding the arm cuff, AC adapter and batteries is warranted to be free from defects in materials and workmanship appearing within 5 years from the date of purchase, when used in accordance with the instructions provided with the monitor. The arm cuff and AC adapter are warranted to be free from defects in materials and workmanship appearing within 1 year from the date of purchase when the monitor is used in accordance with the instructions provided with the monitor. The above warranty extends only to the original retail purchaser.

We will, at our option, replace without charge any monitor, arm cuff or AC adapter covered by the above warranty. Replacement is our only responsibility and your only remedy under the above warranty.

To obtain warranty service contact Customer Service by calling **1-800-634-4350** for the address of the inspection center and the return shipping and handling fee.

Enclose the original printed receipt. Include a letter, with your name, address, phone number, and description of the specific problem. Pack the product carefully to prevent damage in transit. Because of possible loss in transit, we recommend insuring the product with return receipt requested.

THE FOREGOING IS THE SOLE WARRANTY PROVIDED BY OMRON IN CONNECTION WITH THIS PRODUCT, AND OMRON HEREBY DISCLAIMS ANY OTHER WARRANTIES, EXPRESS OR IMPLIED, INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. IMPLIED WARRANTIES AND OTHER TERMS THAT MAY BE IMPOSED BY LAW, IF ANY, ARE LIMITED IN DURATION TO THE PERIOD OF THE ABOVE EXPRESS WARRANTY.

OMRON SHALL NOT BE LIABLE FOR LOSS OF USE OR ANY OTHER SPECIAL, INCIDENTAL, CONSEQUENTIAL OR INDIRECT COSTS, EXPENSES OR DAMAGES.

This warranty provides you with specific legal rights, and you may have other rights that vary by jurisdiction. Because of special local requirements, some of the above limitations and exclusions may not apply to you.

FOR CUSTOMER SERVICE

Visit our web site at:

OmronHealthcare.com

Call toll free:

1-800-634-4350

13. Guidance and Manufacturer's Declaration

OMRON Battery-operated / AC Adapter-operated Blood Pressure Monitor Information for Accompanying Documents in the Scope of IEC60601-1-2:2014

Important information regarding Electromagnetic Compatibility (EMC)

BP7450 manufactured by OMRON HEALTHCARE Co., Ltd. conforms to IEC60601-1-2:2014 Electromagnetic Compatibility (EMC) standard.

Further documentation in accordance with this EMC standard is available at OmronHealthcare.com/emc. Refer to the EMC information for BP7450 on the website.

The **Bluetooth**[®] word mark and logos are registered trademarks owned by the Bluetooth SIG, Inc. and any use of such marks by OMRON HEALTHCARE Co., Ltd. is under license. Other trademarks and trade names are those of their respective owners.

App Store and App Store logo are service marks of Apple Inc., registered in the U.S and other countries.

Google Play and the Google Play logo are trademarks of Google LLC.



OMRON HEALTHCARE Co., Ltd.
53, Kunotsubo, Terado-cho, Muko, Kyoto, 617-0002 JAPAN

Distributed by:
OMRON HEALTHCARE, INC.
1925 West Field Court Lake Forest, IL 60045 U.S.A.

OmronHealthcare.com
© 2019 OMRON HEALTHCARE, INC.

Made in Vietnam